Ci-Ca Dialysate K2 Plus

en Instructions for use

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Fresenius Medical Care AG

Eise-Kröner-Str. 1, 61352 Bad Homburg, GERMANY

☎ +49 6172 609-0

www.freseniusmedicalcare.com/en/product-information

Ci-Ca Dialysate K2 Plus

en

INSTRUCTIONS FOR USE

Please read the following instructions carefully.

GENERAL INFORMATION

General description of the product

Ci-Ca Dialysate K2 Plus is free of bacterial endotoxins and is a sterile, clear and colourless bicarbonate buffered dialysate solution with phosphate for continuous veno-venous haemodialysis (CVVHD) and post-dilution continuous veno-venous haemodiafiltration (postCVVHDF) with citrate for regional anticoagulation (RCA). Each unit consists of a bag with 5,000 mL solution and an injection port and there is a HF connector and a Luer-Lock connector for connection to the tubing system of the blood processing device. The bag system is sterilized using steam. It is a calcium-free dialysis solution.

COMPOSITION

Ci-Ca Dialysate K2 Plus is supplied as a double-chamber bag with a total volume of 5,000 mL. The ready-to-use dialysis solution is obtained by mixing both compartments immediately prior to the application. The small compartment contains 250 mL acidic solution with 20 mmol/L Mg²⁺, 40 mmol/L K⁺, 57 mmol/L Cl⁻, 25 mmol/L H₂PO₄⁻, an H⁺-excess of 2 mmol/L (equivalent to a pH of about 3.1), and 111 mmol/L glucose. The large compartment contains 4,750 mL of a solution with 140 mmol/L Na⁺, 118.84 mmol/L Cl⁻ and 21.16 mmol/L HCO₃⁻.

Each litre of the ready-to-use solution contains:

Sodium chloride	6.604	g/L
Sodium hydrogen carbonate	1.68	g/L
Potassium chloride	0.0559	g/L
Magnesium chloride hexahydrate	0.2033	g/L
Potassium dihydrogen phosphate	0.1701	g/L
Glucose monohydrate	1.1	g/L
Theoretical osmolality	279	mOsm/L
pH	7.0 – 7.6	

Excipients: Water for injections, hydrochloric acid 25%

Ion and glucose concentration in each litre of the ready-to-use solution:

K ⁺	2	mmol/L
H ₂ PO ₄ ⁻	1.25	mmol/L
Na ⁺	133	mmol/L
Mg ²⁺	1	mmol/L
Cl ⁻	115.75	mmol/L
HCO ₃ ⁻	20	mmol/L
Glucose, anhydrous	5.55	mmol/L

INTENDED PURPOSE AND RELATED DEFINITIONS

Intended purpose

Correction of blood electrolytes and acid-base balance in an extracorporeal dialysis treatment.

Medical indication

Acute kidney injury, fluid overload or intoxications requiring CRRT

Intended patient population

Ci-Ca Dialysate K2 Plus is intended for use in adult patients requiring CRRT:

- continuous veno-venous haemodialysis (CVVHD) or
- post-dilution continuous veno-venous haemodiafiltration (postCVVHDF).

Intended user group and intended environment

Ci-Ca Dialysate solutions must only be used by health care professionals who are sufficiently trained in CRRT and in the application of the involved products. The CRRT devices must be operated in an intensive care unit or under similar conditions,

with close medical supervision and continuous monitoring. Ci-Ca Dialysate solutions must be used with CRRT devices that have dedicated pumps for the administration of citrate and calcium and a therapy program supporting RCA, such as Fresenius Medical Care CRRT devices multiFiltrate and multiFiltratePRO.

SIDE EFFECTS

Side effects related to the product and treatment:

Organ class (IMDRF)	Side effects
Nervous system	Seizures, cerebral oedema, demyelination
Blood and lymphatic system	Anaemia*, coagulation disorders (incl. thrombocytopaenia)*, hypovolaemia and hypovolaemic shock (through blood loss, excess ultrafiltration)*, (sustained) hypervolaemia and oedema*, haemolysis*
Immune system	Hypersensitivity/allergic reactions (leading to e.g., dyspnea, hypotension, urticaria, fever and chills, flushes, unconsciousness, abdominal pain)
Vascular system	Bleeding/blood loss* (also catheter-related), air embolism*, thrombosis and thrombo-embolism*
Heart	Cardiac arrhythmia, cardiac arrest
Respiratory system	Bronchospasm, respiratory arrest
Gastro-intestinal system	Vomiting*, abdominal pain*
Metabolism and nutrition	Removal of nutrients*, acidosis (e.g., from citrate accumulation), alkalosis (e.g., from citrate overload), electrolyte imbalance (including hypokalaemia when using low-potassium solutions [K2/K2 Plus], hypophosphataemia and/or hypomagnesaemia, especially when using phosphate-free solutions [K2/K4]), glucose imbalance (e.g., in patients with diabetes mellitus)
Musculo-skeletal system	Cramps/muscle spasms, bone fractures
Infections	Bacterial infection (also: catheter-related) leading to sepsis
Other/generalised disorders	Hypotension*, hypothermia*, headache

*Side effect predominantly related to the treatment mode

For further information on these side effects, see MEDICAL WARNINGS AND PRECAUTIONS as well as INTERACTIONS.

Additional side effects may be specific to other devices or drugs used in the therapy.

Reporting of serious incidents

If any serious incident occurs in relation to the device, including those not listed in this leaflet, the treating physician shall be informed immediately. Within the EU the user and/or patient must report any serious incident that has occurred in relation to the device to the manufacturer according to labelling (M) and the competent authority of the EU Member State in which the user and/or patient is established.

A serious incident can be any incident that directly or indirectly leads to the death of a patient, user or other person; to the temporary or permanent serious deterioration of a patient's, user's or other person's state of health; or a serious public health threat.

MEDICAL WARNINGS AND PRECAUTIONS

(to prevent Side Effects)

- A pre-existing ionized hypocalcaemia may have to be treated before initiating the CRRT procedure to reduce the risk of any clinically relevant hypocalcaemia during the first hours of treatment.
- To avoid severe disturbances, the fluid balance, acid-base status, and serum electrolytes (e.g., Ca²⁺, Na⁺, K⁺, Mg²⁺, inorganic phosphate), as well as blood glucose must be monitored at regular intervals prior to and during treatment. If necessary, the prescription should be adapted, and a severe imbalance must be managed using established medical practices.
- Severe dysnatraemia generally requires a slow serum sodium normalisation, as severe complications could otherwise result, e.g., demyelination or cerebral oedema.
- In patients with reduced citrate metabolism, e.g., in patients with reduced hepatic function, hypoxemia or a disturbed oxygen metabolism, regional citrate anticoagulation can lead to citrate accumulation. Signs include ionized hypocalcaemia, an increased need for calcium substitution, an increased ratio of total ionized calcium (above 2.25), or metabolic acidosis. It may then become necessary to increase the Ci-Ca Dialysate flow, reduce the blood flow, or stop using sodium citrate 4% for anticoagulation and to use an alternative blood anticoagulation agent. Intensified monitoring is recommended.
- Under RCA, loss of filter performance (i.e., clogging) limits removal of uremic toxins and may lead to citrate overload. Signs of citrate overload include metabolic alkalosis, hypernatraemia and hypercalcaemia. In these cases, the extracorporeal circuit should be replaced.
- Patients in a prolonged immobilized position may undergo bone remodelling/demineralisation, ultimately leading to bone fractures. Under RCA, early indications of an ionized hypercalcaemia may be masked by a decrease in the calcium infusion rate. In patients under RCA for longer than 2 weeks or whose calcium infusion rate is progressively decreasing, bone turnover markers should be closely monitored.

Additional medical warnings and precautions may be specific to other devices or drugs used in the therapy.

CONTRAINDICATIONS

Product-specific contraindications

- Hypersensitivity to any constituents of the product
- Hypokalaemia
- Severe hyperphosphataemia

Treatment-related contraindications

- Inability to establish vascular access
- Known severely impaired citrate metabolism (for further information, see MEDICAL WARNINGS AND PRECAUTIONS)

Relative treatment-related contraindications

(predictors for poor treatment outcome – treatment decision on an individual basis)

- Mitochondrial dysfunction potentially leading to impaired citrate metabolism (e.g., paracetamol and metformin intoxications).
- The need to have treatment effects reached more rapidly, e.g., certain intoxications, than the CRRT treatment can accomplish. Intermittent HD may be more appropriate in these cases.
- Terminal disease with no reasonable expectations of recovery.
- Severe dysnatraemias might be better managed with a different anticoagulation method in which it is more straightforward to modify the impact on serum sodium concentration.
- The need to extend the treatment beyond 4 weeks continuously.

Additional contraindications may be specific to other devices or drugs used in the therapy.

Ci-Ca Dialysate K2 Plus

en

PERFORMANCE CHARACTERISTICS AND CLINICAL BENEFITS

Performance characteristics

The dialysis solution must only be used in combination with 4% trisodium citrate infusion (RCA). The reduced bicarbonate and sodium concentration in Ci-Ca Dialysate K2 Plus is concerted to the sodium and bicarbonate intake by trisodium citrate and its metabolites. Ci-Ca Dialysate K2 Plus is calcium-free to sustain the anticoagulatory effect of citrate in the filter. Since it is a calcium-free dialysis solution, the separate infusion of calcium is mandatory.

Ci-Ca Dialysate K2 Plus contains phosphate to restrict the loss of phosphate from the patient's blood during CVHD or postCVHDF therapies.

Clinical benefits

Clinical benefits of CVHD or postCVHDF using RCA in patients suffering from acute kidney injury, fluid overload or intoxications in an intensive care setting may include improved morbidity and improved survival outcomes by gently controlling fluid, by implementing an acid-base and electrolyte balance and by extending therapy time up to continuous application, thereby limiting the rate of changes in the patient (haemodynamic stability with slow fluid status changes and lower risk of cerebral oedema with slow osmotic pressure changes).

INTERACTIONS

- CRRT may reduce the blood concentration of certain medicinal products and nutrients (e.g., those with a low protein binding capacity, with a small distribution volume, and with a molecular weight below the cut-off of the haemofilter). An appropriate revision of the dose of such medicinal products may be required. The removal of important nutrients should be compensated for with adapted (par)enteral nutrition.
- Crystalloids, (par)enteral nutrition and other infusions are commonly given in intensive care medicine. Interactions may be expected with medicinal products that target or show as a side effect the alteration of the blood electrolyte, acid-base content, or fluid status of the patient.
- The cardiac toxicity of cardiac glycosides, in specific digoxin, may be exacerbated upon the correction of hyperkalaemia, hypermagnesaemia, hypocalcaemia, or the development of hyponatraemia or an alkalosis.

Additional interactions may be specific to other devices or drugs used in the therapy.

METHOD OF ADMINISTRATION

During the application, additional devices and solutions are required: A 4% trisodium citrate solution must be infused in predilution at a dose of about 4 mmol citrate per litre of treated blood for CVHD and of about 5 mmol citrate per litre of treated blood for postCVHDF. The citrate dosing should remain within 3 – 5 mmol/L for CVHD and 3 – 5.5 mmol/L for postCVHDF to prevent the occurrence of hypo- or hypernatraemia when Ci-Ca Dialysate K2 Plus is used at the generally recommended ratio to the blood flow (i.e., of 1:3). A calcium solution with a calcium concentration of 50 – 500 mmol/L (preferably 100 mmol/L) must be infused systemically or into the venous bloodline of the CVHD or postCVHDF system immediately prior to the connection with the venous catheter lumen. The infused amount of calcium must be adjusted accordingly to control the systemic ionized calcium (iCa) concentration. A suitable starting dose is normally 1.7 mmol calcium per litre of effluent flow.

The metabolic acid-base status of the patient can be changed by the ratio between the infusion of buffer bases, depending on the blood and citrate solution flow, and the removal of buffer bases with the flow of dialysis solution. It must be noted that an increase of the dialysis solution flow causes a shift towards acidosis. This differs from the application of dialysis solutions (with e.g., 35 mmol/L bicarbonate), which are used with CVHD or postCVHDF in combination with systemic anticoagulation. At a dose of 2,000 mL/h Ci-Ca Dialysate K2 Plus, a balanced metabolic acid-base status is typically achieved in combination with a blood flow of 100 mL/min. This is valid when a sufficiently large high-flux dialyser is used, which does not limit the transport of buffer bases and calcium-citrate complexes.

Dosing

The efficacy of the continuous haemodialysis treatment is generally achieved with an application of 1,500 – 2,500 mL/h dialysis solution in adults, depending on body weight. The daily dose should not exceed 2,000 mL/kg.

Handling Instructions

Before use

Labelling must be clear and the expiry date may not be exceeded.

For convenient handling, the recommended temperature of the solution is +20 °C.

Handling


Opening of the overwrapping

The double-chamber bag containing the dialysis solution must only be taken out of its overwrapping immediately prior to its use. Before opening, the overwrapping must be checked for damage. Bags with damaged overwrapping must be discarded. The peel seam that separates the two compartments must be intact after unpacking; otherwise compliance of the solution specification cannot be guaranteed and the bag must be discarded.

Mixing of both compartments

Both compartments must be mixed immediately prior to the use of the dialysis solution. Do not use unmixed. Mixing of the compartments is mandatory to avoid immediate electrolyte and glucose disturbances during application.

C)

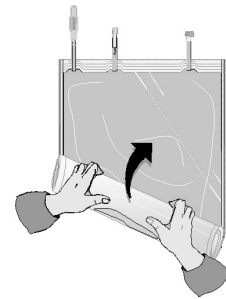... until the peel seam between both compartments has opened along its entire length and the solutions from both compartments are mixed.

After mixing both compartments, a check must be carried out to verify that the peel seam is completely open, that the solution is clear and colourless and that the bag is not leaking. Use only if the solution is clear and the container is not damaged. A leaking bag must be discarded immediately.

Connection of the bag to extracorporeal circuit

Ci-Ca Dialysate K2 Plus is supplied with a yellow HF connector to be used in combination with Fresenius Medical Care CRRT devices multiFiltrate and multiFiltratePRO, and the corresponding Ci-Ca tubing lines containing a compatible yellow HF male connector.

The treating physician is responsible for any use of the Luer-Lock connector on the bag to connect Ci-Ca Dialysate K2 Plus using CRRT devices other than multiFiltrate and multiFiltratePRO.

After mixing both compartments, the bag with the dialysis solution is connected to the dialysate circuit according to the instructions for the chosen device for continuous renal replacement therapy. Contamination of the dialysis solution or any parts in contact with the dialysis solution must be avoided. Any bag of Ci-Ca Dialysate K2 Plus that has been disconnected from the dialysate circuit must be discarded.

Injection port

Use of the injection port is the responsibility of the treating physician and it can be used with established cannula not exceeding 23 Gauge.

Termination


The treatment must be concluded according to the instructions provided with the CRRT device.

WARNINGS AND PRECAUTIONS

Technical Warnings and Precautions


- The RCA protocol must be complied with to limit the occurrence of side effects. Please choose the correct prescription, based on the actual serum blood values for electrolytes and acid base.
- Ci-Ca Dialysate K2 Plus must not be used for intravenous infusion or for direct infusion into the extracorporeal circuit.
- The product is for single use only to minimise the risk of microbial contamination.
- It must be used within 48 hours after mixing to ensure that the chemical properties are maintained.
- Prior to application, Ci-Ca Dialysate K2 Plus must have a temperature of at least +20 °C (also see IFU of CRRT device on temperature of operating room) to support the warming procedure at the CRRT device.
- Depending on clinical requirements, the ready-to-use solution should be warmed immediately before use to 36.5 °C – 38.0 °C. Please adjust the temperature at the CRRT device accordingly.
- The set-up of the Citrate-CVHD or Citrate-postCVHDF must be carefully checked prior to the start of the treatment. Specifically, it is important to ensure the correct set-up of the citrate and calcium infusion. The correct set-up should

A)

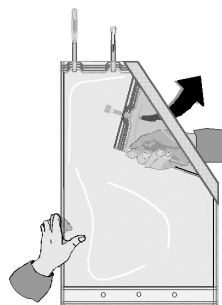Unfold the small compartment.

B)

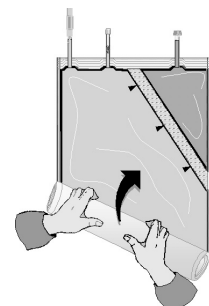Roll up the solution bag starting from the corner opposite the small compartment ...

Ci-Ca Dialysate K2 Plus

en

be confirmed by measuring the citrate-induced drop of the ionized calcium concentration in the extracorporeal circuit no later than 20 – 30 minutes after the start of the treatment. If this drop is missing, the set-up must be re-checked as mixing-up of the citrate and calcium infusion can lead to severe electrolyte imbalance and acid-base disturbance.

- The treating physician is responsible for any substance that is added to the unmixed or mixed solution. The compatibility of the substance must be checked beforehand to avoid undesirable effects such as particle formation. After admixture, the solution must be mixed again.
- Any unused portion of solution is to be discarded.

Additional technical warnings and precautions may be specific to other devices or drugs used in the therapy.

PARTICULAR NOTES ON MATERIALS AND SUBSTANCES

There is no hypersensitivity to the bag system or any of the materials known.

For SVHC information according to Article 33 of Regulation (EC) No. 1907/2006 ("REACH"), please use this page:

www.freseniusmedicalcare.com/en/svhc

SPECIAL PRECAUTIONS FOR STORAGE

Keep away from sharp objects such as knives and scissors.

Do not store below +4 °C or above +25 °C.

Keep dry.

DISPOSAL

Ensure safe disposal of any unused product or waste material in accordance with local regulations.

Materials that have been in contact with blood or other material of human origin may be infectious. Dispose of such materials by taking the necessary precautionary measures and in accordance with local regulations for (potentially) infectious materials.

SYMBOLS USED ON LABELS

	Sterilized using steam
	Catalogue number
	Medical Device
	Manufacturer
	Use-by date
	Single sterile barrier system with protective packaging outside
	Patient information website
	Units
	Temperature limit

	Batch code
	Unique Device Identifier
	Date of manufacture
	Caution
	Do not re-use
	Consult instructions for use

DATE OF REVISION OF TEXT

2022-05

A digital version of this document is available on the following website:

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