

# SENSA® i - maximum comfort and flexibility

LIKAMED Product Portfolio



**FRESENIUS  
MEDICAL CARE**

# Introduction

## LiKAMED

LiKAMED shares more than 25 years of successful partnership and expertise with Fresenius Medical Care in the field of medical patient chairs and bed chairs and is one of the leading providers of dialysis chairs in Europe and abroad. Both companies share a common vision of offering innovative products to improve the patient's quality of life. Products from LiKAMED fulfill the highest standards in terms of patient comfort, ergonomics, functionality and user safety.

Made by LiKAMED means: all development, design and manufacturing is done at the factory in Germany. All chair types are thoroughly tested and additionally subject to voluntary testing and certification by an independent TÜV institute.

The medical patient chairs and bed chairs from LiKAMED have a strong international position in our NephroCare and customer clinics and are fully integrated into our Fresenius Medical Care logistics

and technical service processes. LiKAMED products, with the associated technical service, allow Fresenius Medical Care to offer a better range of products and services from one single source, thereby strengthening our relationship with our customers, providing them with solutions and creating value.

### LiKAMED company certifications

LiKAMED products are manufactured in Germany, tested according to applicable norms and, additionally, carrying TÜV product test mark. Through years of experience and rigorous quality controls, we stand together for highest quality „Made in Germany“.

LiKAMED is certified

- QM-certification according to ISO 9001
- QM-certification according to EN ISO 13485
- Certified Environmental Management System according to EN ISO 14001



### Our partner:

#### LiKAMED GmbH

For the past 40 years, LiKAMED has been highly successful in the medical sector. The formula for success: design and quality. The idea: Using this formula and many years of experience to develop innovative products for practical – and convenient – use. The claim: All products are



developed, tested, manufactured and assembled by LiKAMED in-house specialists and are certified according to the applicable ISO management systems. The promise: Made by LiKAMED – Made in Germany – Focused on Design and Quality.

# SENSA® i – high flexibility and function

The new SENSА® i offers very high standards for configuration, safety, comfort and ease of use. It continues the tradition of LiKAMED quality in a new design marked by complete flexibility. In addition to its high quality, close attention has been paid to the comfort of patients and users, as shown by the wide selection of configurable equipment options.

All Sensa® i dialysis chairs from Fresenius Medical Care are equipped with equipotential bonding cable and full foam armrests as well as Trendelenburg position and 7.5 cm casters with individual locking are included as standard options.

- up to five silent drives for enjoyable and stepless adjustment of the chair
- flexible armrests for every application
- wide range of comfortable upholsteries and colours
- multiple types of covers
- footrests can be adjusted manually or electrically
- central locking system for quickly moving or locking the chair
- multiple types of casters with different functions and floor clearances
- equipotential bonding
- safe working load 250 kg



## Configure your chair exactly as required

The modular design allows you to assemble the chair exactly as you need it. You configure your chair precisely to your individual requirements with different function groups such as casters, covers, armrests, upholsteries and much more.



Stepless motorized adjustment from seat-, comfort- and bed position to shock- or Trendelenburg position

# Flexibility

## Adjustment and positioning options

Set the chair to any desired position. Depending on the adjustment options and the installed options, the SENSE® i chair comes with three, four or five drives.



### Flexibility with seat position, comfort and bed position

Stepless adjustable positioning of back part, seat part and foot part with three motors.



### Additional option with height adjustment

Particularly for older patients it is extremely important to be able to get on or off the chair easily, which is why the additional motor for stepless adjustment of the height is so important. This also allows the chair to be moved into a specific working position, which is a great benefit for the nursing personnel.

## Shock or Trendelenburg position

You can configure your chair to select whether it moves the patient into the shock position or the Trendelenburg position in case of emergency. Trendelenburg position is the standard, shock position can be chosen instead as an option.



*shock position*



*Trendelenburg position*



## Ease of operation

All adjustments are easily made quickly and precisely with the intuitively designed hand switch. The symbols on the buttons are three-dimensional for operation by visually impaired people. You receive the suitable hand switch for the specific model with only the functions that you have selected.

### Comfort buttons

One of the most popular functions of LiKAMED chairs: after setting the optimal seat position, the patient can change the position quickly and easily with only one pair of buttons without changing the angles of the seat part, back part and foot part to one another.

### Special buttons

The coloured markings ensure that the most important functions are easily visible. Quick adjustment to the shock or Trendelenburg position is particularly important in case of emergency.

### Magnetic key

The hand switch can be locked and unlocked with the magnetic key.

### Foot controller

Foot controller in an emergency, the optional foot controller can be double-clicked to move the chair quickly and easily into the shock or the Trendelenburg position, depending on your configuration. At the same time, nursing personnel have their hands free for immediate action in case of emergency.





# Comfort

The SENSE® i offers a wide range of different upholsteries including the perfect selection for your application. The wide range of optional equipment, such as spring core, soft upholstery or individually


adjustable air chambers with dynamic adjustment options, make it easy for you to find the best solutions for your application.

## Upholsteries



### Air chambers

The hardness of the seat and back parts can be adjusted for additional comfort during long treatment sessions.

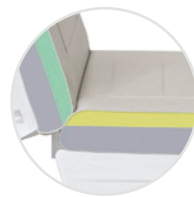

### Soft

The soft upholstery has a particularly soft-top layer.



### Visco foam

The visco foam is particularly yielding and adjusts itself to the contour of the patient by the body heat, thereby enabling comfort while seated and lying down for many hours.



### Universal

The two-layer sandwich foam is ideal for universal applications with many different patients.



### Spring core

The inner spring core has a particularly long life and a slightly firmer comfortable seat while seated and lying down.




### Robust

The PU foam is particularly suitable for heavier patients, allowing them to be seated or lie down in comfort.

## Individualised colour design

All cushion covers can be configured in single colour or two-tone in the seven standard colours at no additional cost. A total of 37 custom colours is also available for additional colour configuration.

The Fresenius Medical Care standard is stone and sapphire.



Stone



Saffron



Raspberry



Azure



Kiwi



Agave



Sapphire

## Pulsating seat – Perfect for treatments that take place over several hours


The firmness of this upholstery can be adjusted to suit the patient with internal air chambers.



The air chambers in the back part and the seat part can be adjusted quickly and quietly with the separate hand switch to make it as comfortable as possible for the patient. The ability to adjust the air chambers for patients of different weights is also a great advantage when treating a range of different patients.



The switchable pulsating movements of the back part relieve the spine and relax the musculature. These small and continuous adjustments provide enjoyable patient comfort and prevent cramped sitting position due to static sitting.

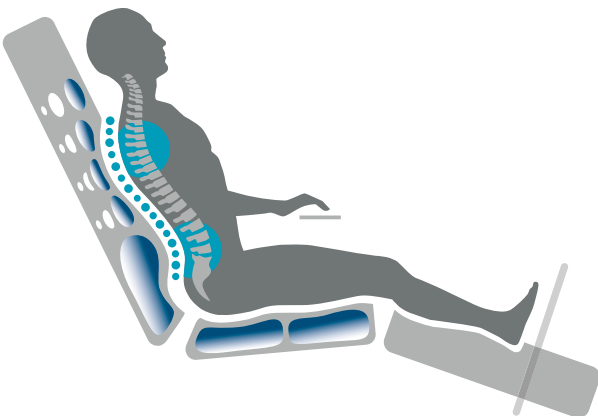

- Separate adjustment possibilities of the back part for individual comfort
- Pulsating adjustment of the back part for relaxation and relief of the spine
- Large air chambers in the seating area to correct the angle of the pelvis and reduce the pressure of the intervertebral discs
- Adjustable sitting hardness according to individual patient weight

# Comfort

## Relax pillow and neck roll

A relax pillow or a neck roll can be used for additional comfort in the chair.




**Relax pillow**  
*For a relaxed head position.*



**Neck roll**  
*For relaxation in the neck region.*

## Armrests

For every specialty a good solution: armrests that are rotatable and can be folded upwards, armrests that are tiltable or armrests with height adjustment.

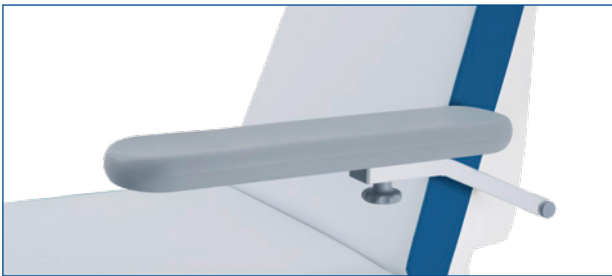

**Rotatable and folding upwards armrests**  
*The folding armrests make it easier for patients to take their position. The rotatable function provides increased comfort while seated or lying down.*



**Additional tiltable armrests**  
*Increased comfort even during long treatment sessions due to extended adjustment options. For excellent positioning for dialysis treatment.*



**Additional height adjustable and tiltable armrests**  
*Excellent adjustment of position in all seated and lying positions. Greater comfort and flexibility for particularly large or small patients.*



## Armrest pads

A wide range of different pads for armrests makes it easy to make the right choice for your purpose.

The full foam armrests pads are comfortable and particularly easy to clean and therefore the Fresenius Medical Care standard. The range of foamed armrest pads with imitation leather covers is particularly comfortable. Arms are comfortably supported by the


slightly concave shape and the integrated holder gives the hand switch its specific place and it is easily accessible.

The imitation leather cover on the foam armrest pads enables you to select different colours to design the chair to your requirements and to match the interior design of its surroundings.



### **Full foam armrest pads**

*The grey, full-foamed armrest pads are soft, noted for their long-lived and ease of cleaning.*

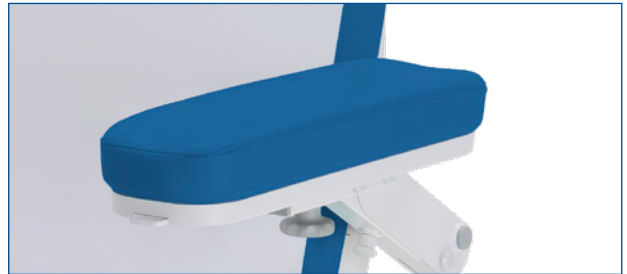

### **Foam armrest pads**

*The straight armrest pads with a holder for the hand switch are particularly comfortable with the soft padding material.*



### **Low wedge-shaped foam armrest pads**

*Very good if you want an armrest position with a slight inclination downwards for greater comfort.*



### **High wedge-shaped foam armrest pads**

*This type improves comfort during treatment with a strong inclination downwards, even with not tiltable armrests.*

# Comfort

## Footrest

The footrest provides a really relaxed seat position. A footrest is recommended for specific applications, such as dialysis, in order to give patients the option of pressing their feet against the footrest to relieve a cramp.



## Footrest



### **Without footrest**

*Treatment chairs without footrest are especially space saving.*



### **Manually adjustable footrest**

*Greater stability and comfort for patients, particularly during extended treatment times. Nursing personnel can adjust the footrest with the setting wheel on the side.*



### **Automatic footrest length compensation**

*When the foot part is adjusted, the length compensation moves the footrest slightly outwards when it is raised and slightly inwards when it is lowered. This ensures that legs are never confined.*



### **Electrically adjustable footrest with length compensation**

*The patient can easily adjust the position of the footrest with the hand switch. This can save time for personnel during the treatment period. The automatic length compensation is also integrated with the electrical footrest.*


## Casters

Whether the chair is moved around often or is always in the same place, we have the right casters for every situation.



### **Casters Ø 7.5 cm with individual lock**

*A reduced access height simplifies entry and exit for the patient and assistance by the nursing personnel. The individual lock on all four casters keeps it stationary. If the chair is not moved around much, the caster diameter of 7.5 cm is a good solution.*

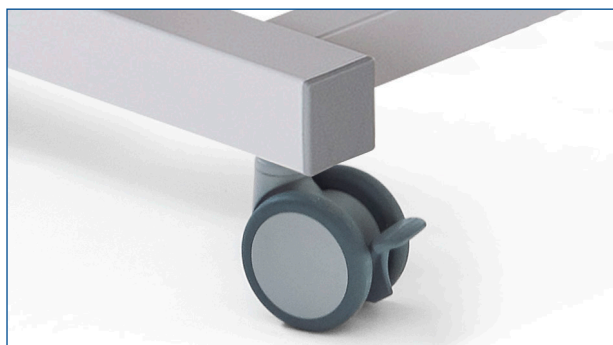

### **Casters Ø 10 cm with individual lock**

*Easy to move combined with individual lock if the chair is not moved very often.*



### **Casters Ø 10 cm with individual lock and direction locking**

*Easy to move with good shock absorption with 10 cm diameter casters and individual locking. If the chair is moved over long distances or through corridors, the direction can be locked. It locks the direction of a caster and prevents the chair from moving to the side when it is transporting.*

## Casters with central lock



### **Casters Ø 10 cm with central lock**

*Easy to move around and comfortable with good shock absorption with 10 cm diameter casters. The integrated central locking greatly simplifies setup and adjustment of the chair, particularly in a confined space. The central lock casters includes the additional COMPLETE underframe cover.*



### **Casters Ø 12.5 cm with central lock and direction locking**

*If the chair is moved frequently, a central lock for quick release and braking in combination with the 12.5 cm diameter caster is recommended. The direction locking fixes the direction of a caster and helps to prevent the chair from moving to the side when it is transporting. This makes it easy to maintain the direction of the chair, particularly over longer distances. This model also provides greater floor clearance for accessories and includes the additional COMPLETE underframe cover. The higher seating position makes the work of the nursing personnel easier and more comfortable.*

## Push handles



### **Push handle**

*The push handle makes the chair easier to move and offers the patient more comfort without hands on the head part.*



### **Rotatable push handles**

*When folded out, the push handles make it easier to move the chair, even with heavy patients. They can also be used when folded in and save space in confined areas.*



## Battery

Rechargeable battery for emergency situations to be able to operate in case of power cuts. Use the chair adjustment options without a power supply.

# Hygiene and cleaning

## Covers

The different covers make cleaning simple and easy and save time for other purposes. You can fit covers over the seat part, the foot part and the bottom frame. The back part cover is part of standard equipment.

The cover components are practical and also look attractive as design elements.

## Upholstery protection

The robust and easily cleaned transparent protective covers protect your chair even under heavy use with shoes. We also recommend a Jersey cover as a colour accent and protection for your chair. The hand switch holder on the armrests remains accessible with all protective covers.



### **Back part cover**

*Smooth, easily cleaned surface. All Sensa® i chairs come with the elegant back part cover with a practical holder for bags and accessories or for the power cord when moving the chair.*



### **Seat part cover**

*The optional seat part covers offer quick and simple cleaning and a visual design variant.*



### **Foot part rear cover**

*The cover protects the bottom of the foot part and greatly simplifies cleaning.*



### **Column cover**

*This elegant and easily cleaned cover encloses the vertical column or the lift column with the height adjustment.*



### **Underframe cover LIGHT**

*The column area cover is very easy to clean and fits all chair models with individually locking casters.*



### **Underframe cover COMPLETE**

*This underframe cover is an excellent option for fastest and easiest cleaning. It fits all types of casters and is standard equipment on all chair models with central locking.*



## Transparent protective covers

Easily cleanable and available for foot part, footrest and all armrest pads.



## Antimicrobial upholstery

The special antimicrobial imitation leather is available for all upholsteries and can reduce the viability or infectiousness of microorganisms, thus offering greater safety for patient and health care personnel.



## Paper roll holder

Hygiene and cleanliness for every patient by changing the paper cover. The permanently installed roll holder keeps the paper roll ready for use and positions the paper cover.



# Accessories

## SENSA® i accessories



### **LED reading lamp**

*Flexible positioned light for patient comfort in a low-light environment.*




### **Extension adapter for LED reading lamp**

*For very tall patients. Extends the space for working around the armrest.*



### **Lumbar support**

*The individually adjustable cushion promotes the natural curve of the spine in a seated position and may prevent backache.*

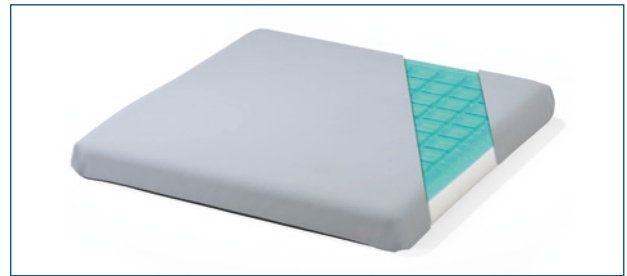

### **Anti-decubitus gel seat cushion**

*For use with extended treatment sessions and for a special seated experience.*



### **Hand switch holder with velcro**

*Made of imitation leather and suitable for full-foam armrest pads. Makes it easier for the patient to grip the switch.*




### **Infusion rod holder**

*For improved mobility by space-saving mounting swiveling infusion rod on the right and left of the chair.*



### **Infusion rod incl. bottle holder**

*The height-adjustable accessory for the swiveling infusion rod holder.*

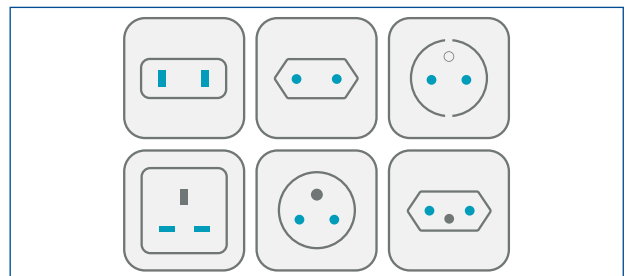

### **Power plug types**

*Country-specific power plug.*

# Equipment and options – configuration options

## Overview of variations

Flexibility		
Adjustability	back, seat and foot part motorized adjustable	●
	height motorized adjustable	○
Emergency positions	shock position	○
	Trendelenburg position	●
	foot controller for shock or Trendelenburg position	○
Patient comfort		
Upholsteries	universal 2-layer foam upholstery	●
	adjustable air chambers upholstery	○
	soft upholstery	○
	visco foam upholstery	○
	spring core upholstery	○
	robust upholstery	○
Headrest	relax pillow	○
	neck roll	○
Armrests	rotatable and folding upwards armrests	●
	additional tiltable armrests	○
	tiltable and height adjustable armrests	○
Armrest pads	full foam armrest pads	●
	foam armrest pads	○
	low wedge-shaped foam armrest pads	○
	high wedge-shaped foam armrest pads	○
Footrests	manually adjustable footrest	○
	manually adjustable footrest with automatic length compensation	○
	electrically adjustable footrest with automatic length compensation	○
Mobility		
Casters	casters Ø 7.5 cm with individual lock	●
	casters Ø 10 cm with individual lock	○
	casters Ø 10 with individual lock and direction locking	○
	casters Ø 10 cm with central lock	○
	casters Ø 12.5 cm with central lock and direction locking	○
Push handles	push handle	○
	rotatable push handles	○
Electric power independence	battery	○

● = standard    ○ = option

# Equipment and options – configuration options

## Overview of variations

Hygiene and cleaning		
Covers	seat part cover	<input type="radio"/>
	foot part bottom cover	<input type="radio"/>
	underframe cover LIGHT	<input type="radio"/>
	underframe cover COMPLETE (standard with central locking)	<input type="radio"/>
	column cover	<input type="radio"/>
Upholstery protection	antimicrobial upholstery	<input type="radio"/>
	transparent, cleanable protective covers for armrest pads	<input type="radio"/>
	transparent, cleanable protective cover for foot part	<input type="radio"/>
	transparent, cleanable protective cover for footrest	<input type="radio"/>
	jersey stretch cover for upholstery	<input type="radio"/>
	jersey stretch cover for relax pillow or neck roll	<input type="radio"/>
	jersey stretch cover for armrest pads	<input type="radio"/>
	paper roll holder	<input type="radio"/>
Accessories		
	LED reading lamp	<input type="radio"/>
	extension adapter for LED reading lamp	<input type="radio"/>
	adjustable lumbar support	<input type="radio"/>
	anti-decubitus gel seat cushion	<input type="radio"/>
	hand switch holder (with velcro for full foam armrest pads)	<input type="radio"/>
	infusion rod holder	<input type="radio"/>
	infusion rod incl. bottle holder	<input type="radio"/>
	equipotential bonding cable	<input type="radio"/>
	country-specific power plug types	<input type="radio"/>

● = standard    ○ = option



## How to configure your chair

The LiKAMED product configurator guides you to the right configuration step by step.

[https://fmc.configurator.likamed.de/en\\_US/configurator/](https://fmc.configurator.likamed.de/en_US/configurator/)



# Technical Specifications

Technical data	SENSA® i
Armrests length x width x height	
full foam armrest pads	60.5 cm x 13.5 cm x 4.5 cm
foam armrest pads	52.5 cm x 15.5 cm x 8.5 cm
low wedge-shaped foam armrest pads	52.5 cm x 15.5 cm x 6.5 to 9.5 cm
high wedge-shaped foam armrest pads	52.5 cm x 16.5 cm x 7.5 to 12.5 cm
Upholstery length x width x height	219 cm x 60 cm x 9 cm
Total length (max.) x total width (max.)	229 cm x 87 cm
Floor clearance	63 mm (casters Ø 7.5 cm) / 93 mm (casters Ø 10.0 cm) / 132 mm (casters Ø 12.5 cm)
Seating height	
with height lift	53 cm to 75 cm (casters Ø 7.5 cm) / 56 cm to 78 cm (casters Ø 10.0 cm) / 60 cm to 82 cm (casters Ø 12.5 cm)
without height lift	57 cm (casters Ø 7.5 cm) / 60 cm (casters Ø 10.0 cm) / 64 cm (casters Ø 12.5 cm)
Height adjustment	22 cm (optional)
Angle of inclination	
back/seat/foot part	0° to 72° (shock position model), -12° to 72° (Trendelenburg position model) / 0° to 25° / -25° to 27°
back/seat/foot part in shock position	0° / 25° / 27°
back/seat/foot part in Trendelenburg position	-12° / 25° / 27°
Equipotential bonding terminal	standard
Power rating, primary	100-240 V AC, 50/ 60 Hz
Power rating, secondary	29 V
Fuses	internal
Protection type hand control unit / controller	IP X6/IP X4
Protection class	II
Duty cycle	2 min. ON / 18 min. OFF or 10% with max. 5 switching cycles per minute
Consumption in standby	< 2 W
Electrical safety, EMC	DIN EN 60601-1, 60601-1-2
Spiral power cable	2.0 m to 3.5 m
Permitted air pressure	800 to 1,060 hPa
Ambient temperature	10° to 40° C
Emission sound pressure level	< 65 dB(A)
Medical device class	I
Applied parts type	B (CF limit values in compliance)
Safe working load	250 kg
Covers and metal frame colour	light grey RAL 7047
Weight (depending on configuration)	85 to 100 kg
Guarantee	2 years

# Upholstery colours

## Standard colours



Stone



Saffron



Raspberry



Azure



Kiwi




Agave



Sapphire

## Special colours

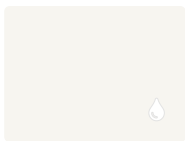

White



Cream



Lemon



Orange



Red



Pastel Blue



Lime



Chestnut



Light Grey



Chrome



Anthracite



Black

◆ Also available as antimicrobial upholstery.





**Additional colours for applications, neck roll, relax pillow, lumbar support,  
all foam armrests**



White



Cream



Pearl



Stone



Sand



Saffron



Corn



Lemonade



Lemon



Mandarin



Orange



Rosewood



Heath



Cyclame



Raspberry



Red



Cherry



Mint



Turquoise



Agave



Kiwi



Lime



Amazonas



Petrol



Violet



Orchid



Melange



Terra



Nut



Chestnut



Light Blue



Azure



Bleu



Baltic



Sapphire



Atoll



Royal



Marine



Chrome



Slate



Anthracite



Black



Pastel Blue

**Features**

The imitation leather of the upholstery is highly durable. The upholstery is tearproof and resistant to dirt, scratches and disinfectants.

# Fresenius Medical Care packages for Sensa® i

## The Sensa® i includes always:

- 3 motors for adjustment of back-, foot- and seat part
- Trendelenburg position
- grey full foam armrests
- individually locking casters 7.5 cm
- equipotential bonding cable

## Available packages



### Sensa® i Package Basic

- Manual footrest



### Sensa® i Package Value

- height motorized adjustment
- armrest tiltable & height adjustable
- 10 cm casters with central locking
- manual footrest with autom. length adjustment



### Sensa® i Package Premium

- with height motorized adjustment
- armrest tiltable & height adjustable
- 10 cm casters with central locking
- electrical footrest with autom. length adjustment
- relax pillow with Fresenius Medical Care logo
- visco foam upholstery
- foot part rear cover

# DT 1 and ST 1 – Side tables

## Features

- Single-handed, spring supported height adjustment
- Removable tabletop with finger-hollow
- Tilttable tray (only model DT 1)
- Detachable plastic storage box
- Holder, e.g. for bags
- Strong metal frame construction
- Plastic-coated subframe with excellent durability, protection against corrosion and abrasion

## Options



Holder for monitor etc., Ø 16 mm  
incl. clamping screw (DT 1 and ST 1)\*



Holder for monitor etc., Ø 16 mm  
incl. clamping screw (only ST 1)\*



Lid for storage box, dark grey



Container/bottle holder

- Flat large tray, not removable, fixed to the frame (only ST 1)



Side table DT 1, tiltable tray incl. storage box



Side table ST 1, fixed tray incl. storage box

Model	DT 1 and ST 1
Height	84–98 cm (adjustable)
Area required	42 × 66 cm
Castors	Ø 5 cm
Storage box:	1
Capacity	6 l
Dimensions	appr. 17 × 27 × 15 cm
Colour	Light Grey RAL 7047
Tray:	2
Size	40 × 53 cm (removable)
Small, trapezoidal	40/30 × 18 cm (fixed)
Colour	Light Grey RAL 7047
Colour of frame	Light Grey RAL 7047
Maximum load	25 kg
Weight	17 kg

\*Only available for new orders.

# Article numbers

## Upholstery colour cover / sides

### Sensa® i Dialysis Chairs

F00008903	<b>LIKAMED – Sensa® i</b> (with EU-power plug line, equipotential bonding cable, 3 motors for adjustment of back-, foot- and seat part, Trendelenburg position, grey full foam armrests, individually locking casters 7,5 cm diameter, upholstery surface colour in stone, upholstery side colour in sapphire, <b>NO</b> pillow/neck roll, <b>NO</b> foot rest)	Stone/Sapphire
-----------	--	----------------

### Sensa® i Packages – FME Standard packages (for new orders only)

F00009303	<b>LIKAMED Sensa® i package Basic</b> (with manual footrest F00008912)
F00009304	<b>LIKAMED Sensa® i package Value</b> (with high motorized adj. (F00008904), armrest tiltable & height adjustable (H260261), 10 cm casters with central locking (F00008917), manual footrest incl. autom. length adjustment (F00008913))
F00009305	<b>LIKAMED Sensa® i package Premium</b> (with high motorized adj. (F00008904), armrest tiltable & height adjustable (H260261), 10 cm casters with central locking (F00008917), electrical footrest incl. autom. length adjustment (F00008915), relax pillow with FME logo (F00008907), visco foam upholstery (5073801), foot part rear cover (F00009057))

### Sensa® i Options: Adjustability and Emergency Positions (for new orders only)

F00008904	<b>LIKAMED Sensa® i height motorized adjustable</b>
F00008905	<b>LIKAMED Sensa® i with shock position</b> (instead of Trendelenburg position)
5073771	<b>LIKAMED Sensa® i foot controller for emergency</b> (Shock or Trendelenburg) <b>position</b>

### Sensa® i Upholstery & Head-, Arm-, Footrest Options (for new orders only)

F00008906	<b>LIKAMED Sensa® i adjustable air chambers upholstery</b>
H260471	<b>LIKAMED Sensa® i with special soft seat part</b> (Warranty period for this upholstery is 6 months)
5073801	<b>LIKAMED Sensa® i visco foam upholstery</b>
5073821	<b>LIKAMED Sensa® i springs core at seat part</b>
5073811	<b>LIKAMED Sensa® i upholstery robust seating area</b>
5073931	<b>LIKAMED chair antibacterial surface upholstery</b> (Only available in conjunction with these colours: Stone, Saffron, Raspberry, Azure, Kiwi, Creme, Lemon, Orange, Red, Sapphire, Lime, Chestnut, Black, Melange and Anthracite)
F00008907	<b>LIKAMED Sensa® i relax pillow with FME logo</b> (In upholstery side colour)
F00008908	<b>LIKAMED Sensa® i neck roll with FME logo</b> (In upholstery side colour)
5074671	<b>LIKAMED chair pillow / neck roll without FME logo</b>
5073951	<b>LIKAMED chair lumbar support</b> (In upholstery cover colour)
F00008432	<b>LIKAMED chair gel seat cushion grey</b>

**Sensa® i Upholstery & Head-, Arm-, Footrest Options** (for new orders only)

F00008909	<b>LIKAMED Sensa® i two foam armrests</b> (In upholstery colour, cover and sides)
F00008910	<b>LIKAMED Sensa® i two low wedge-shaped foam armrests</b> (In upholstery colour, cover and sides)
F00008911	<b>LIKAMED Sensa® i two wedge-shaped foam armrests</b> (In upholstery colour, cover and sides)
H260351	<b>LIKAMED Sensa® i two armrests tiltable</b>
H260261	<b>LIKAMED Sensa® i two armrests tiltable and height adjustable</b>
F00008912	<b>LIKAMED Sensa® i manual footrest</b> (Inlay in upholstery sides colour)
F00008913	<b>LIKAMED Sensa® i manual footrest with length compensation</b> (Inlay in upholstery sides colour)
F00008915	<b>LIKAMED Sensa® i electrical footrest incl. length compensation</b> (without foot part rear cover) (Inlay in upholstery sides colour)

**Sensa® i Mobility & Frame Cover Options** (for new orders only)

5073991	<b>LIKAMED Sensa® i individual locking casters 10 cm diameter</b>
F00008916	<b>LIKAMED Sensa® i individual locking casters 10.0 cm diameter and direction locking</b>
F00008917	<b>LIKAMED Sensa® i casters 10.0 cm diameter with central locking – CA</b>
F00008918	<b>LIKAMED Sensa® i casters 12.5 cm diameter with CA and direction locking</b>
H260271	<b>LIKAMED Sensa® i push handles, fixed</b>
5073761	<b>LIKAMED Sensa® i push handles, rotatable</b>
H260201	<b>LIKAMED Sensa® i rechargeable battery</b>
F00008919	<b>LIKAMED Sensa® i seat part covers</b>
F00008920	<b>LIKAMED Sensa® i underframe cover LIGHT</b>
F00008921	<b>LIKAMED Sensa® i underframe cover COMPLETE</b> (standard with central locking casters)
F00008922	<b>LIKAMED Sensa® i column cover</b>
F00009057	<b>LIKAMED Sensa® i foot part rear cover</b>

**Further Optional Accessories** (for new orders only)

H260221	<b>LIKAMED Sensa® i paper roll holder</b>
H260211	<b>LIKAMED Sensa® i LED reading lamp</b>
H260971	<b>LIKAMED extension adapter for LED lamp</b>
H260231	<b>LIKAMED Sensa® i infusion rod holder</b>
H260241	<b>LIKAMED Sensa® i infusion rod with bottle holder</b> (Only available in conjunction with Sensa® i infusion rod holder)
H260861	<b>LIKAMED hand switch holder Velcro grey for full foam armrest</b>
F00008841	<b>LIKAMED hand switch holder Velcro grey <u>long</u> for foam armrest</b>

# Article numbers

**Sensa® i Colour Options** (for new orders only) Please note: special colours have an additional cost and may cause longer delivery time for the order

Article No.	Upholstery cover	Article No.	Upholstery cover
5074691	<b>Stone</b> (the standard colour in 2020 and after)	F00008533	<b>Light grey</b> (the standard colour 2019, 2020+ special)
5074701	<b>Saffron</b> (standard colour)	5074711	<b>Raspberry</b> (standard colour)
5074721	<b>Azure</b> (standard colour)	5074781	<b>Sapphire</b> (standard colour)
5074611	<b>Agave</b> (standard colour)	5074731	<b>Kiwi</b> (standard colour)
5074681	<b>Pastel blue</b> (special colour)	5074151	<b>White</b> (special colour)
5074741	<b>Cream</b> (special colour)	5074751	<b>Lemon</b> (special colour)
5074761	<b>Orange</b> (special colour)	5074771	<b>Red</b> (special colour)
5074791	<b>Lime</b> (special colour)	5074801	<b>Chestnut</b> (special colour)
5074811	<b>Black</b> (special colour)	F00008442	<b>Chrome</b> (special colour)
F00008472	<b>Anthracite</b> (special colour)		
Article No.	Upholstery sides	Article No.	Upholstery sides
F00008928	<b>Stone</b> (the standard in 2020 and after)	F00008935	<b>Light grey</b> (the standard colour 2019, 2020+ special)
F00008929	<b>Saffron</b> (standard colour)	F00008930	<b>Raspberry</b> (standard colour)
F00008932	<b>Azure</b> (standard colour)	F00008931	<b>Sapphire</b> (standard colour)
F00008933	<b>Agave</b> (standard colour)	F00008934	<b>Kiwi</b> (standard colour)
F00008980	<b>Pastel blue</b> (special colour)	F00008981	<b>White</b> (special colour)
F00008982	<b>Cream</b> (special colour)	F00008983	<b>Lemon</b> (special colour)
F00008984	<b>Orange</b> (special colour)	F00008985	<b>Red</b> (special colour)
F00008986	<b>Lime</b> (special colour)	F00008987	<b>Chestnut</b> (special colour)
F00008988	<b>Black</b> (special colour)	F00008989	<b>Chrome</b> (special colour)
F00008990	<b>Anthracite</b> (special colour)		
<b>LIKAMED Sensa® i upholstery sides in another special colour</b> _____			
Please find possible further special colours in product brochure and note colour name here.			
F00008924	<b>LIKAMED chair lumbar support in same colour as upholstery sides</b>		
F00008925	<b>LIKAMED chair relax pillow/neck roll in same colour as upholstery cover</b>		
F00008926	<b>LIKAMED chair relax pillow/neck roll strap in same colour as upholstery sides</b>		
F00008927	<b>LIKAMED chair foot rest inlay in same colour as upholstery cover</b>		



**Likamed Power Plug Lines**

Please note: If not otherwise indicated dialysis chairs will be delivered with EU plug line

H261011	<b>LIKAMED power plug line Switzerland</b>
H260281	<b>LIKAMED power plug line GB</b>
H260251	<b>LIKAMED power plug line USA MCLII</b>
H261021	<b>LIKAMED power plug line Israel</b>
H260361	<b>LIKAMED power plug line Brazil</b>
	<b>LIKAMED power plug line other</b> (Please indicate)

**Sensa® i Upholstery Cover Options** (for new orders only)

5074101	<b>LIKAMED chair jersey cover for upholstery</b>
H260991	<b>LIKAMED chair jersey stretch cover for armrest</b> (One piece)
5074141	<b>LIKAMED chair jersey cover for wedge armrest</b> (One piece)
5074131	<b>LIKAMED chair jersey cover for low wedge armrest</b> (One piece)
5074111	<b>LIKAMED chair jersey cover for neck roll</b>
5074121	<b>LIKAMED chair jersey cover for relax pillow</b>
H260311	<b>LIKAMED chair transparent cover for full foam armrest</b> (One piece)
5074091	<b>LIKAMED chair transparent cover wedge armrest</b> (One piece)
5074081	<b>LIKAMED chair transparent cover low wedge foam armrest</b> (One piece)
F00008975	<b>LIKAMED Sensa® i transparent cover foot part</b>
F00008923	<b>LIKAMED Sensa® i transparent cover footrest</b>

**Likamed Side Tables & Options** (for new orders only)

6296881	<b>LIKAMED – ST 1 – Side table, fixed tray</b> (incl. storage box)
6296891	<b>LIKAMED – DT 1 – Side table, tiltable top</b> (incl. storage box)
H260491	<b>LIKAMED DT1/ST1 bottle holder, integrated in the fixed small tray</b>
5074881	<b>LIKAMED ST1 flat large tray not removable, fixed to the frame</b>
5074891	<b>LIKAMED ST1 holder for monitor etc., fixed to the large tray-frame</b>
5074901	<b>LIKAMED DT1/ST1 holder for monitor etc., fixed to small tray-frame</b>
5074911	<b>LIKAMED DT1/ST1 lid for storage box, dark grey colour</b>

**Manufacturer:**

LiKAMED GmbH  
Raußmühlstraße 7  
75031 Eppingen  
Germany  
Phone: +49 7262 9189-0  
Fax: +49 7262 9189-900  
E-Mail: [info@likamed.de](mailto:info@likamed.de)  
[www.likamed.de](http://www.likamed.de)



Sensa is protected as trademark in selected countries

Some equipment is not available in all configurations. The information is up to date at the time of printing. Fresenius Medical Care Deutschland GmbH accepts no liability for discrepancies of colour and shape from the pictures. No responsibility is accepted for errors and misprints. Fresenius Medical Care Deutschland GmbH reserves the right to make changes to publications and products at any time.

Where the masculine form or the generic masculine is used in this brochure for the sake of legibility, this specifically includes the feminine and other gender identities.



**FRESENIUS  
MEDICAL CARE**

Head office: Fresenius Medical Care Deutschland GmbH · 61346 Bad Homburg v. d. H. · Germany  
Phone: +49 (0) 6172-609-0 · Fax: +49 (0) 6172-609-2191  
[www.FreseniusMedicalCare.com](http://www.FreseniusMedicalCare.com)